

# ZEPTOLink IOL Positioning System™

Designed to enhance the precision and efficiency of cataract surgery, the ZEPTOLink IOL Positioning System™ is the first and only device that creates an instantaneous capsulotomy for precise, 360-degree IOL overlap and improved outcomes.



Centered capsulotomies



360° overlap



Improved outcomes


## ZEPTOLink™ POWER CONSOLE:

**Seamless integration** with any phacoemulsification system

**Small ergonomic design** and user-friendly touchscreen display

**Efficient** system setup, prime, and operation

**Provides complete surgeon control** via phaco foot pedal

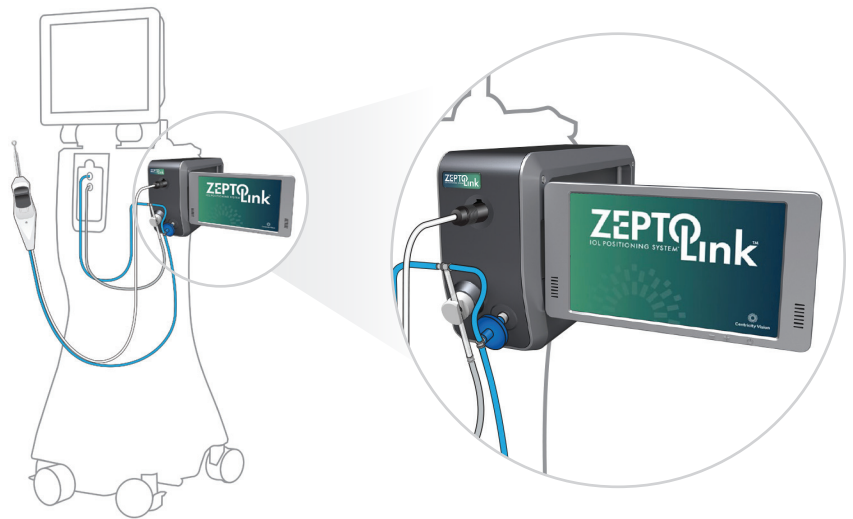

## HANDPIECES:

Nitinol ring design delivers **consistent energy and temperature uniformity** for strong and precise 360° capsulotomies

Suction provided by the phacoemulsification system allows for **precise, controlled, consistent, and centered capsulotomies**

**ZEPTOLink™ Handpiece and ZEPTOLink Plus™ Handpiece offer 2 capsulotomy sizes** to better support the surgical needs of each patient, ideal for premium IOL and complex cases

## ZEPTOLink™ HANDPIECE



**ZEPTOLink Plus™ HANDPIECE**

For more information, visit [zeptozone.com](http://zeptozone.com) or call **1-833-339-3786**.

# TECHNICAL SPECIFICATIONS

## ZEPTOLink IOL Positioning System™

### ZEPTOLink™ POWER CONSOLE

#### Parameter

#### Specification

Power input	100-240 volts AC, 50-60 Hz, 1.66 amps
Energy delivery	1.13 joule
Vacuum requirement of OEM ultrasonic phacoemulsification system (UPS)	Minimum suction of 535 mm Hg
Case Dimensions	12.5 in (L) x 4.8 in (W) x 8.0 in (H)
Weight	8.0 lb

### ZEPTOLink™ HANDPIECE

#### Parameter

#### Specification


Average capsulotomy diameter	5.1 mm +/- 0.08 mm
Cutting element	Nitinol ring, Ø 4.4 mm
Silicone suction cup	Ø 5.9 mm, height 1.17 mm

### ZEPTOLink Plus™ HANDPIECE

#### Parameter

#### Specification

Average capsulotomy diameter	5.4 mm +/- 0.17 mm
Cutting element	Nitinol ring, Ø 4.7 mm
Silicone suction cup	Ø 6.2 mm, height 1.17 mm

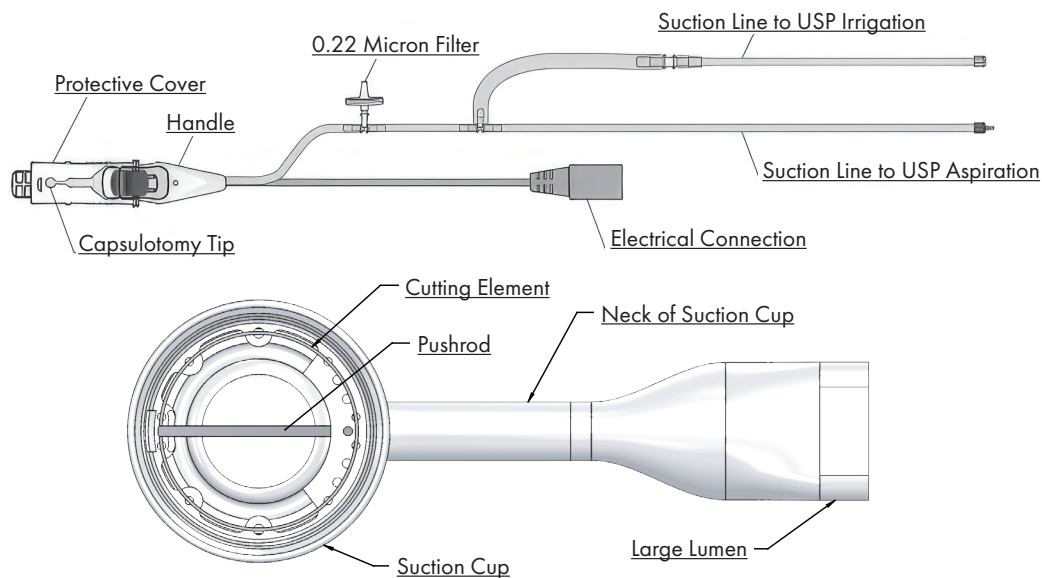

For more information, visit [zeptozone.com](http://zeptozone.com)  
or call **1-833-339-3786**.